

LOOK OF A STROLLER. TECHNOLOGY OF A WHEELCHAIR.





Portable Design

Every element of Flip was designed to encourage parents and their family to get out and go. With a simple "flip" and fold like a stroller, you're on your way. Little Wave Flip transforms into a small package that can fit in the trunk of the tiniest car. Now there's room for all the other essentials needed for your adventure.










Find yours at www.kimobility.com/splashzone

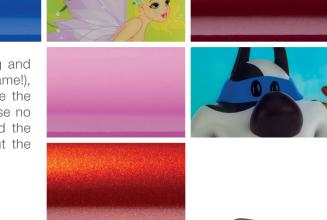





Don't struggle with growth, embrace it. Adjusting the seat height, depth and width has never been so easy to change. More built in adjustment and fewer parts needed for growth than any other chair.



Contracture Footrest

Allows for extension, flexion, or lower leg rotation for optimal angle of knees, ankles and windswept deformity. Configure it bilateral or center mount with the same parts.



Peds Pro ELR

Patented adjustable pivot location to promote proper knee placement and an adjustable calf pad for personalized fit.

SPECIFICATIONS



Multi-Angle Footrest

Infinite adjustment in multiple axes to accommodate difficult positioning requirements such as inversion, eversion, plantar, and dorsiflexion. Optimizes support, pressure distribution, and postural alignment.



Help shield your child from the sun's harmful rays or when the sky is dark and gray. The unique Flip-Up visor channels water off to the side or extends downward for extra protection from the elements.



Attendant Wheel Lock

A hands-free, simple to use stroller like wheel lock option that makes a parent's life just a little easier.

Transport Weight:	21 lbs.*	9.53 kg*
Seat Width:	10-18 in.	25-46 cm
Seat Depth:	12-20 in.	30-51 cm
Seat Height:	13.5-20 in.	34-51 cm
Tubing Size:	1 in.	2.54 cm
Weight Capacity:	165 lbs.	75 kg
WC-19 Transit:	Yes	Yes

^{* 14&}quot;x14" Frame with Attendant Wheel Locks. Armrests, Rear Wheels and Footrests Removed.

