

Durable medical equipment **product guide**



Innovative products and support



Wheelchairs

When choosing a wheelchair, there are many factors to consider including Arm Style, Leg Style, Weight Capacity and Seat Dimension. We've issued the charts in a way that makes it easier to select the right one.

Arm Style: Many users require the armrests in order to safely perform transfers in and out of a wheelchair. Individuals who can stand also use the armrests for support when getting in and out of the wheelchair. Properly fitted armrests when used help in maintaining a correct and functional seated posture. Arm Style choices include:

- Desk-length swing back and removable
- Desk-length swing-back and height adjustable
- Desk-length height-adjustable swing-back and removable
- Full-length removable
- Full-length removable and height adjustable
- Full-length swing-back and removable
- Full-length height-adjustable swing-back and removable

Leg Style: The main purpose of these leg rests is to redistribute the pressure on the user and increase comfort. There are two options for leg style:

- Swing-Away
- Elevating

Weight Capacity: Wheelchairs come with varying weight limits. Basic wheelchairs have weight limits of 250–350 lbs, lightweight or pediatric wheelchairs have weight limits of 200–250 lbs, and extra wide and bariatric wheelchairs have weight limits of 500–850 lbs.

Seat Dimensions: Seat widths typically range from 16 to 20 inches, with an average width of about 18 inches. Take seat depth into consideration as it will make a difference in overall comfort and usage. Seat depth is the distance from the back of the seat to the front edge and should be measured when the user is sitting straight up.




Wheelchair Sizing Guide

Determine the best fit for your wheelchair in six easy steps.

The most important step to consider when choosing a wheelchair is a proper fit. It can make the difference between recovering well and causing an injury.

Improper fit can cause:

- Posture issues
- Breathing problems
- Sores
- Discomfort
- Pelvic issues
- Hip problems

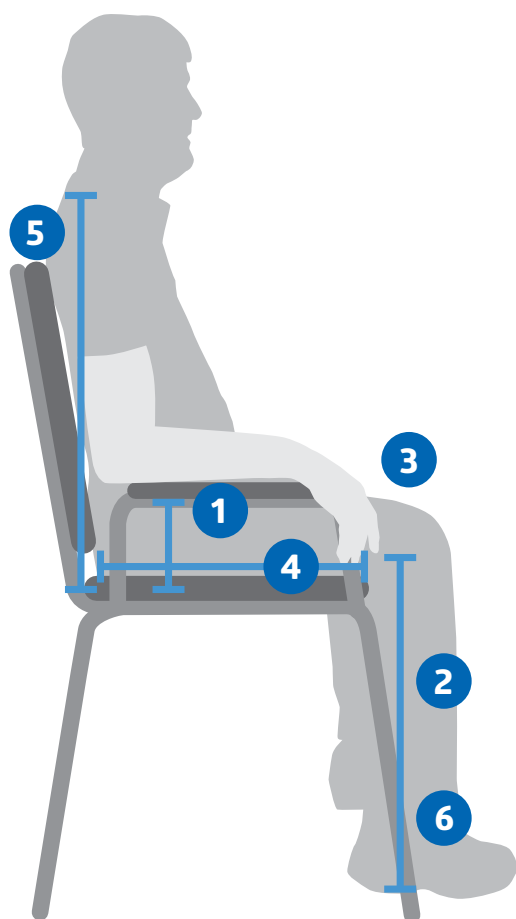

1 Determine Arm Type and Height

- **Full-length arm:** Gives support to pivot transfers when the patient pushes off to stand.
- **Desk-length arm:** Makes it easy to move around tables and desks.
- **Height:** Measure from the elbow to the seat of the chair while holding arms up with elbows bent at a 90-degree angle. Use adjustable arms when possible.

2 Determine Footrest Style

- **Elevated leg rest:** Raises legs
- **Articulation leg rest:** Extends longer for taller patients.
- **Length:** Measure from back of the knee to the heel of the foot.

3 Determine Wheelchair Weight and Weight Limit

- **Weight capacity:** Choose wheelchair weight based on the patient's weight.
- **Upper body strength:** Choose lighter wheelchairs for weaker patients.

4 Determine Seat Size

- **Seat width:** Measure from hip to hip in a straight line and then add two inches.
- **Seat depth:** Measure from the back of the hip to the back of the knee while seated and then subtract one inch.

5 Determine Back Height

- **Standard back height:** Measure from collarbone to the seat while seated.
- **Taller back height:** Require additional upper trunk support to be installed.

6 Determine Floor to Seat Height

- **Movement:** Check if the patient uses feet to propel or move.
- **Height:** Measure the distance between the back of the knee and the heel.

K1 Guardian Wheelchairs

- Durable hammertone frame; available in nylon or vinyl upholstery
- Hemi-height adjustable
- Accommodates 2" seat extender

Nylon

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length Swing Back	Swing-Away	K1166N22S	Black	300 lb (136.1 kg)	16" x 16"	K0001 OR K0002 + E0951	N	1 ea
		K1186N22S	Black	300 lb (136.1 kg)	18" x 16"	K0001 OR K0002 + E0951	N	1 ea
		K1206N22S	Black	300 lb (136.1 kg)	20" x 16"	K0001 OR K0002 + E0951	N	1 ea
	Elevating	K1166N22E	Black	300 lb (136.1 kg)	16" x 16"	K0001 OR K0002 + E0990	N	1 ea
		K1186N22E	Black	300 lb (136.1 kg)	18" x 16"	K0001 OR K0002 + E0990	N	1 ea
		K1206N22E	Black	300 lb (136.1 kg)	20" x 16"	K0001 OR K0002 + E0990	N	1 ea
Full-Length Removable	Swing-Away	K1186N13S	Black	300 lb (136.1 kg)	18" x 16"	K0001 OR K0002 + E0951	N	1 ea
	Elevating	K1186N13E	Black	300 lb (136.1 kg)	18" x 16"	K0001 OR K0002 + E0990	N	1 ea

Vinyl

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length Swing Back	Swing-Away	K1166V22S	Black	300 lb (136.1 kg)	16" x 16"	K0006	N	1 ea
		K1186V22S	Black	300 lb (136.1 kg)	18" x 16"	K0006	N	1 ea
		K1206V22S	Black	300 lb (136.1 kg)	20" x 16"	K0006 + E2201	N	1 ea
	Elevating	K1166V22E	Black	300 lb (136.1 kg)	16" x 16"	K0006 + K0195	N	1 ea
		K1186V22E	Black	300 lb (136.1 kg)	18" x 16"	K0006 + K0195	N	1 ea
		K1206V22E	Black	300 lb (136.1 kg)	20" x 16"	K0006 + E2201 + K0195	N	1 ea
Full-Length Removable	Swing-Away	K1186V11S	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea
	Elevating	K1186V11E	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea



Swing-Away
Arm Rest

Accommodates
2" seat depth
extenders

K1166V22S

K2 Guardian Wheelchairs

- Durable Hammertone frame
- Ships in hemi height; includes seat belt, anti-tippers and brake extenders

Nylon

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Full-Length Removable and Height Adjustable	Elevating	K2166N42E1	Black	300 lb (136.1 kg)	16" x 16"	K0002 + E0973 + K0195	N	1 ea
		K2186N42E1	Black	300 lb (136.1 kg)	18" x 16"	K0002 + E0973 + K0195	N	1 ea
		K2206N42E1	Black	300 lb (136.1 kg)	20" x 16"	K0002 + E0973 + K0195	N	1 ea

New Guardian Extra-Wide Chairs

- Durable Hammertone frame and wipeable vinyl upholstery
- Heavy-duty frame with 500lb weight capacity
- Hemi-height adjustable

Vinyl

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Removable Desk-Length Arms	Swing-Away	KX208V21S	Black	500 lb (226.8 kg)	20" x 18"	-	N	1 ea
		KX228V21S	Black	500 lb (226.8 kg)	22" x 18"	-	N	1 ea
		KX248V21S	Black	500 lb (226.8 kg)	24" x 18"	-	N	1 ea
		KX268V21S	Black	500 lb (226.8 kg)	26" x 18"	-	N	1 ea
	Elevating	KX208V21E	Black	500 lb (226.8 kg)	20" x 18"	-	N	1 ea
		KX228V21E	Black	500 lb (226.8 kg)	22" x 18"	-	N	1 ea
		KX248V21E	Black	500 lb (226.8 kg)	24" x 18"	-	N	1 ea
		KX268V21E	Black	500 lb (226.8 kg)	26" x 18"	-	N	1 ea

Accommodates 2" seat extenders



K2166N42E1

Includes anti-tippers



MDS806550PLUS

The K3 wheelchair comes standard with the accessories your customers are looking for.

New K3 Guardian Wheelchairs

- Durable Hammertone frame and easy-to-clean nylon upholstery
- Swing-back arms are also removable for easy transfers and storage
- Hemi-height adjustable

Height-adjustable swing-back arms for comfort and easy side transfers

Seat belt prevents sliding or falling out of the wheelchair

Choose from elevating or swing away foot rests

2" (5.1cm) seat depth extender for customer comfort

Break extenders allow the user to independently lock the breaks



K3186N44E1



Anti-tippers for added safety



Heel loops for support and to prevent sliding (swing-away footrest models only)

New K3 Guardian Wheelchairs

- Durable Hammertone frame and easy-to-clean nylon upholstery
- Swing-back arms are also removable for easy transfers and storage
- Hemi-height adjustable

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length Swing-Back and Removable	Swing-Away	K3166N24S	Black	300 lb (136.1 kg)	16" x 16"	K0006	N	1 ea
		K3186N24S	Black	300 lb (136.1 kg)	18" x 16"	K0006	N	1 ea
		K3206N24S	Black	300 lb (136.1 kg)	20" x 16"	K0006 + E2201	N	1 ea
	Elevating	K3166N24E	Black	300 lb (136.1 kg)	16" x 16"	K0006 + K0195	N	1 ea
		K3186N24E	Black	300 lb (136.1 kg)	18" x 16"	K0006 + K0195	N	1 ea
		K3206N24E	Black	300 lb (136.1 kg)	20" x 16"	K0006 + K0195 + E2201	N	1 ea
Desk-Length, Height- Adjustable, Swing-Back and Removable	Swing-Away	K3166N34S	Black	300 lb (136.1 kg)	16" x 16"	K0006 + E0973	N	1 ea
		K3186N34S	Black	300 lb (136.1 kg)	18" x 16"	K0006 + E0973	N	1 ea
		K3206N34S	Black	300 lb (136.1 kg)	20" x 16"	K0006 + E2201 + E0973	N	1 ea
	Elevating	K3166N34E	Black	300 lb (136.1 kg)	16" x 16"	K0006 + K0195 + E0973	N	1 ea
		K3186N34E	Black	300 lb (136.1 kg)	18" x 16"	K0006 + K0195 + E0973	N	1 ea
		K3206N34E	Black	300 lb (136.1 kg)	20" x 16"	K0006 + K0195 + E2201 + E0973	N	1 ea
Full-Length, Swing Back and Removable	Swing-Away	K3166N14S	Black	300 lb (136.1 kg)	16" x 16"	K0006	N	1 ea
		K3186N14S	Black	300 lb (136.1 kg)	18" x 16"	K0006	N	1 ea
		K3206N14S	Black	300 lb (136.1 kg)	20" x 16"	K0006 + E2201	N	1 ea
	Elevating	K3166N14E	Black	300 lb (136.1 kg)	16" x 16"	K0006 + K0195	N	1 ea
		K3186N14E	Black	300 lb (136.1 kg)	18" x 16"	K0006 + K0195	N	1 ea
		K3206N14E	Black	300 lb (136.1 kg)	20" x 16"	K0006 + K0195 + E2201	N	1 ea
Full-Length, Height- Adjustable, Swing Back and Removable	Swing-Away	K3166N44S	Black	300 lb (136.1 kg)	16" x 16"	K0006 + E0973	N	1 ea
		K3186N44S	Black	300 lb (136.1 kg)	18" x 16"	K0006 + E0973	N	1 ea
		K3206N44S	Black	300 lb (136.1 kg)	20" x 16"	K0006 + E2201 + E0973	N	1 ea
	Elevating	K3166N44E	Black	300 lb (136.1 kg)	16" x 16"	K0006 + K0195 + E0973	N	1 ea
		K3186N44E	Black	300 lb (136.1 kg)	18" x 16"	K0006 + K0195 + E0973	N	1 ea
		K3206N44E	Black	300 lb (136.1 kg)	20" x 16"	K0006 + K0195 + E2201 + E0973	N	1 ea
		K3166N44E1*	Black	300 lb (136.1 kg)	16" x 16"	K0006 + K0195	N	1 ea
		K3186N44E1*	Black	300 lb (136.1 kg)	18" x 16"	K0006 + K0195	N	1 ea
		K3206N44E1*	Black	300 lb (136.1 kg)	20" x 16"	K0006 + K0195 + E2201	N	1 ea

*Includes anti-tippers, seat extension, brake extension and seat belt.



K4 Lightweight Wheelchairs

- Durable, lightweight TIG-welded frame folds down for easy transport

Nylon

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Swing Back	Swing-Away	MDS806500N	Black	300 lb (136.1 kg)	16" x 16"	K0004	N	1 ea
		MDS806500	Black	300 lb (136.1 kg)	18" x 16"	K0004	N	1 ea
		MDS806560	Black	350 lb (158.8 kg)	20" x 16"	K0004+E2201	N	1 ea
	Elevating	MDS806570	Black	350 lb (158.8 kg)	22" x 16"	K0004 + E2201	N	1 ea
		MDS806550N	Black	300 lb (136.1 kg)	16" x 16"	K0004 + K0195 or E0990	N	1 ea
		MDS806550	Black	300 lb (136.1 kg)	18" x 16"	K0004	N	1 ea
		MDS806565	Black	350 lb (158.8 kg)	20" x 16"	K0004+E2201+K0195	N	1 ea
Full-Length Removable	MDS806575	Black	350 lb (158.8 kg)	22" x 16"	K0004 + E2201 + K0195	N	1 ea	
	Swing-Away	MDS806500FLA	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea
	Elevating	MDS806550FLA	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea

K4 Basic Wheelchairs

- Durable Grey frame
- Hemi height adjustable

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
		MDS806500E	Silver	300 lb (136.1 kg)	18" x 16"	K0004	N	1 ea
		MDS806500NE18	Silver	300 lb (136.1 kg)	16" x 18"	K0004	N	1 ea
		MDS806560E	Silver	300 lb (136.1 kg)	20" x 18"	K0006 + E2201	N	1 ea
	Elevating	MDS806550NE	Silver	300 lb (136.1 kg)	16" x 16"	K0004 + K0195	N	1 ea
		MDS806550E	Silver	300 lb (136.1 kg)	18" x 16"	K0004+K0195	N	1 ea
		MDS806500E18	Silver	300 lb (136.1 kg)	18" x 18"	-		1 ea
		MDS806565E	Silver	300 lb. (136.1 kg)	20" x 18"	K0003 or K0006 + E2201 + K0195 or E0990	N	1 ea



MDS806560



MDS806500E

K4 Plus Wheelchairs

- Durable Hammertone frame
- Includes anti-tippers

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Swing-Back and Height Adjustable	Swing-Away	MDS806500PLUS	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea
		MDS806560PLUS	Black	300 lb (136.1 kg)	20" x 18"	K0004 + E2201 + E0971	N	1 ea
	Elevating	MDS806550PLUS	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea

New Guardian K4 Wheelchairs

- Durable Hammertone frame and easy to clean nylon upholstery
- Hemi-height adjustable
- Quick-release rear wheels

Nylon

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Swing Back and Removable	Swing-Away	K4166N24S	Black	300 lb (136.1 kg)	16" x 16"	-	N	1 ea
		K4186N24S	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea
		K4206N24S	Black	300 lb (136.1 kg)	20" x 16"	-	N	1 ea
	Elevating	K4166N24E	Black	300 lb (136.1 kg)	16" x 16"	-	N	1 ea
		K4206N24E	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea
		K4186N24E	Black	300 lb (136.1 kg)	20" x 16"	-	N	1 ea
Desk-Length, Height-Adjustable, Swing Back and Removable	Swing-Away	K4186N34S	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea
	Elevating	K4186N34E	Black	300 lb (136.1 kg)	18" x 16"	-	N	1 ea



K4186N24S

Guardian Recliners

- Infinite position reclining up to 180°
- Durable hammertone frame

Nylon

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Removable	Elevating	KR168N21E	Black	300 lb (136.1 kg)	16" x 18"	K0006 + E0990 + E1226	N	1 ea
		KR188N21E	Black	300 lb (136.1 kg)	18" x 18"	K0006 + E0990 + E1226	N	1 ea
		KR208N21E	Black	300 lb (136.1 kg)	20" x 18"	K0006 + E0990 + E1226	N	1 ea
Full-Length Removable	Elevating	KR188N11E	Black	300 lb (136.1 kg)	18" x 18"	K0006 + E0990 + E1226	N	1 ea

Reclining Wheelchairs

- Infinite position reclining from 90°–140°
- Durable carbon steel frame
- Hemi-height adjustable

Vinyl

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Removable	Elevating	MDS808350	Chrome	300 lb (136.1 kg)	16" x 17"	K0003 + E0966 + E1226 + K0195 + E0971	N	1 ea
		MDS808450	Chrome	300 lb (136.1 kg)	18" x 17"	K0003 + E0966 + E1226 + K0195 + E0971	N	1 ea
		MDS808550	Chrome	300 lb (136.1 kg)	20" x 17"	K0003 + E0966 + E1226 + K0195 + E0971 + E2201	N	1 ea
		MDS808650	Chrome	300 lb (136.1 kg)	22" x 17"	K0007 + E1226 + E2201 + E0971 + E0966 + K0195 or E0990	N	1 ea

Hybrid Wheelchairs

- Combination wheelchair and transport chair
- Quick release large wheel

Nylon

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Removable	Swing-Away	MDS806250NH2	Blue	300 lb. (136.1 kg)	16" x 16"	K0004 or E1038	N	1 ea
		MDS806250H2	Blue	300 lb (136.1 kg)	18" x 16"	K0001 or K0002	N	1 ea
	Elevating	MDS806300H2	Blue	300 lb (136.1 kg)	18" x 16"	K0004 or E1038 + K0195 or E0990	N	1 ea



Excel Wheelchairs

- Durable carbon steel frame
- Hemi-height adjustable

Vinyl

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Size	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Removable	Swing-Away	MDS806250N	Chrome	300 lb (136.1 kg)	16" x 16"	K0001 or K0002	N	1 ea
		MDS806250D	Chrome	300 lb (136.1 kg)	18" x 16"	K0001 or K0002 + E0990 or K0195	N	1 ea
		MDS806400	Chrome	300 lb (136.1 kg)	20" x 16"	K0001 + E2201	N	1 ea
	Elevating	MDS806300N	Chrome	300 lb (136.1 kg)	16" x 16"	K0001 or K0002 + K0195	N	1 ea
		MDS806300D	Chrome	300 lb (136.1 kg)	18" x 16"	K0001 or K0002	N	1 ea
Full-Length, Permanent	Swing-Away	MDS806450	Chrome	300 lb (136.1 kg)	20" x 16"	K0001 + E2201 + K0195	N	1 ea
	Elevating	MDS806150D	Chrome	300 lb (136.1 kg)	18" x 16"	K0001 or K0002	N	1 ea
		MDS806200N	Chrome	300 lb (136.1 kg)	16" x 16"	K0002 + K0195 or E0990	N	1 ea
Full-Length, Removable	Elevating	MDS806100D	Chrome	300 lb (136.1 kg)	18" x 16"	K0001 or K0002	N	1 ea
		MDS806400FLA	Chrome	300 lb (136.1 kg)	20" x 16"	-	N	1 ea
	Elevating	MDS806300DFLA	Chrome	300 lb (136.1 kg)	18" x 16"	K0001 or K0002	N	1 ea

Excel Extra-Wide Wheelchairs

- Durable carbon steel frame
- Hemi-height adjustable

Vinyl

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Removable	Swing-Away	MDS806700	Chrome	500 lb (226.8 kg)	20" x 18"	K0007	N	1 ea
		MDS806800	Chrome	500 lb (226.8 kg)	22" x 18"	K0007	N	1 ea
		MDS806900	Chrome	500 lb (226.8 kg)	24" x 18"	K0007+E2202	N	1 ea
	Elevating	MDS806750	Chrome	500 lb (226.8 kg)	20" x 18"	K0007 + K0195 or E0990	N	1 ea
		MDS806850	Chrome	500 lb (226.8 kg)	22" x 18"	K0007 + K0195 or E0990	N	1 ea
Full-Length, Removable	Swing-Away	MDS806950	Chrome	500 lb (226.8 kg)	24" x 18"	K0007 + E2202 + K0195	N	1 ea
		MDS806700FLA	Chrome	500 lb (226.8 kg)	20" x 18"	K0007	N	1 ea
		MDS806800FLA	Chrome	500 lb (226.8 kg)	22" x 18"	K0007	N	1 ea
	Elevating	MDS806900FLA	Chrome	500 lb (226.8 kg)	24" x 18"	K0007 + E2202	N	1 ea
		MDS806750FLA	Chrome	500 lb (226.8 kg)	20" x 18"	K0007 + K0195 or E0990	N	1 ea
		MDS806850FLA	Chrome	500 lb (226.8 kg)	22" x 18"	K0007 + K0195 or E0990	N	1 ea
		MDS806950FLA	Chrome	500 lb (226.8 kg)	24" x 18"	K0007 + E2202 + K0195	N	1 ea



MDS806300D



MDS806900

Shuttle Bariatric Wheelchairs

- Durable carbon steel frame
- Hemi-height adjustable

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-length, removable	Elevating	MDS809650	Black	850 lb (385.5 kg)	24" x 20"	K0007 + K0195 or E0990 + E2202 + E2203	N	1 ea
		MDS809750	Black	850 lb (385.5 kg)	26" x 20"	K0007 + E2202 + E2203 + K0195 or E0990	N	1 ea
		MDS809850	Black	850 lb (385.5 kg)	28" x 20"	K0007 + K0108 + K0195 or E0990 + E2203	N	1 ea
		MDS809950	Black	850 lb (385.5 kg)	30" x 20"	K0007 + K0108 + K0195 or E0990 + E2203	N	1 ea

Pediatric Wheelchairs

- Adjustable telescoping handles
- Hemi-height adjustable

Vinyl

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Swing-Back	Elevating	MDS806140PD	Yellow	250 lb (113.4 kg)	14" x 12"	E1236 + E0971 + E0990 or K0195	N	1 ea

Nylon

Arm Style	Leg Style	Item No.	Frame Color	Wt. Capacity	Seat Dimensions	HCPCS	Retail Pkg.	Pkg.
Desk-Length, Swing-Back	Swing-Away	KPD2N22S	Yellow	250 lb (113.4 kg)	12" x 12"	E1236	N	1 ea
		KPD4N22S	Yellow	250 lb (113.4 kg)	14" x 12"	E1236	N	1 ea
		KPD2N22S1	Blue	250 lb (113.4 kg)	12" x 12"	E1236	N	1 ea
		KPD4N22S1	Blue	250 lb (113.4 kg)	14" x 12"	E1236	N	1 ea



MDS809750



MDS806140PD



KPD2N22S



KPD4N22S1

Wheelchair Accessories

Item No.	Description	Pkg.
MDS85190	Oxygen Tank Holder and IV Pole for Wheelchair	1 ea
WCA6100	Oxygen Tank Holder for Adult Wheelchair, Holds D and E Tanks	1 ea
MDS85189	Rear Anti-Tip Device for Excel 2000 Wheelchair	1 pr
STMSM22	Anti-Rollback Device	1 ea
MDS85189H	Rear Anti-Tip Device for K1 Basic Wheelchair	1 pr
MDS85181U	Oxygen Tank Holder for K2 Basic, Excel 2000 and Extra-Wide Wheelchairs, Holds Standard Size Tanks	1 ea
MDS85183	Telescoping IV Pole for Excel Extra-Wide Wheelchair	1 ea
MDS85189WK4	Anti-Tip Rear Wheel for K4 and Recliner	1 pr
WCABELT	Wheelchair Safety Belt	1 ea
MDS85157PR	Black Handle Extension Wheel Lock 5/8" x 8"	1 pr
MDS85197A	Adjustable Anti-Theft Device for Wheelchair	1 ea
MDSK41RMDRIVE	Single-Arm Drive attachment for K4 Basic	1 ea
WCA806991	Tinnerman Nut for Wheelchair	1 pr
WCA6200	Adjustable Cup Holder	1 ea
MDS85196	Anti-Fold Device	1 ea
MDS85195	Super Hemi Height Kit for Wheelchairs with Standard Wheels	1 set
MDS851965SH	Anti-Fold Device for Shuttle Extra-Wide Wheelchair	1 ea
MDS85197SH	Overhead Antitheft Device for Shuttle Extra-Wide Wheelchair	1 ea
MDS85180	5-in-1 IV/O2/Anti-Theft Device	1 ea
STMSM22W	Anti-Rollback Device	1 ea
MDS85189SH	Rear Anti-Tipper for Shuttle Extra-Wide Wheelchair	1 pr
MDS85189K	Rear Anti-Tipper for Kidz Pediatric Wheelchair	1 pr
WCAA0060	K1/K2 Nylon 2" Seat Depth Extender	1 ea
WCAA0061	K1/K2 Nylon 2" Seat Depth Extender	1 ea
WCAA0070	K1 & K2 Anti-Tipper	1 pr
WCAA0069	K3 Anti-Tipper	1 pr
WCAA0066	K3 16" Seat Extender	1 ea
WCAA0067	K3 18" Seat Extender	1 ea
WCAA0068	K3 20" Seat Extender	1 ea



MDS85183